
# DT-100 DIGITAL THERMOMETER

## EN Instructions

The photographs and other images of the product in this manual and on the packaging are accurate as possible, but do not claim to be a perfect reproduction of the product.


### THERMOMETER (FIG. 1)
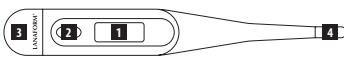
### DISPLAY (FIG. 2)

### INSTALLING THE BATTERY (FIG. 3)
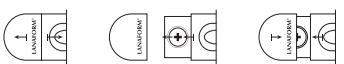
## INTRODUCTION

Thank you for choosing the LANAFORM Digital Thermometer DT-100. Combining accuracy and ease of use, this thermometer allows your baby's rectal, oral or under-arm temperature to be taken. The thermometer lends itself well to measurements being taken at home. With its ultra-flexible tip, the Digital Thermometer DT-100 is comfortable and gentle for your baby. Quick to use, it takes your temperature in just 10 seconds.

## TITLE AND DESCRIPTION OF DISPLAY SYMBOLS ON THE SHEET DESCRIBING THE PRODUCT'S FEATURES

Symbol definition	Identification	Information provided by the manufacturer
Manufacturer's details		Hangzhou Universal Electronic Co., LTD 38 Yangjiatang, Sandun 310030 Westlake District, Hangzhou China
Batch number		Production batch number
BF device		
Instructions for use		Follow the instructions on how to use the device
Protective rating for liquids and solids		Coding system showing the degrees of protection obtained by the electrical equipment being used, preventing access to dangerous parts and penetration by foreign solid objects and water. Provides additional information about this form of protection.
Please note: see instructions for use		Please note: consult the enclosed documents. Please read these instructions before using the device. Please also follow these instructions carefully when using the device.
Notified body		
European representative		Name: Prolinx GmbH Add: Brehmstr. 56, 40239 Dusseldorf Germany Tel: 0049-2131 40519680
Electrical and electronic waste (WEEE)		The device should not be disposed of with other household waste at the end of its serviceable life.

## PLEASE READ ALL THE INSTRUCTIONS BEFORE USING YOUR THERMOMETER, ESPECIALLY THIS BASIC SAFETY INFORMATION

- When using this product, please make sure you follow all the instructions listed below. Any deviation from these instructions may cause injury or affect the accuracy of the readings.
- Contact a doctor if you have any concerns about the accuracy of the readings. Making any diagnosis or starting any treatment without a doctor's advice may pose a risk.
- Always keep the thermometer away from heat sources such as stoves, radiators, etc.
- Do not use this device in a room in which aerosols (sprays) are used or in a room in which oxygen is administered.
- Do not place the thermometer in direct sunlight, or expose it to high temperatures, humidity or dirt.
- Do not immerse the device in water.
- Do not drop the device and do not use it if it is damaged.
- This device must be used at an ambient temperature between 10°C/50°F and 35°C/95°F.
- Do not expose the thermometer to electric shocks.
- Do not expose this thermometer to extremes of temperature above 55°C/131°F or below -25°C/-13°F.
- Do not use the device if the relative humidity is above 80% or below 25%.
- Do not use the device in the vicinity of large electric-magnetic fields, such as those of cordless or portable telephones, etc.
- Please wait at least 30 minutes after taking a shower, exercising, eating or drinking before taking any temperature.
- It is recommended that the thermometer should be kept at room temperature for around 10 minutes before use.
- If there is a problem with your device, please contact your reseller.
- Do not try to repair the device yourself.
- Do not take the device apart as a whole. This device is not intended for use by people, including children, with limited physical, sensory or mental abilities, or by people lacking experience or knowledge, unless they are supervised or have been instructed beforehand on how to use the equipment by someone responsible for their safety. Children should be supervised to make sure they do not play with the device.
- Only use the thermometer for its intended purpose, as described in the instructions.
- The thermometer is pre-set when it leaves the factory. You do not therefore have to calibrate the device when you first use it.
- Body temperature is generally lower in the morning and gradually rises throughout the day. Please therefore take your body temperature at the same time during the day and under the same conditions of use.
- Consult your doctor if you are taking any medicines. These might increase your body temperature.
- The device's performance may be adversely affected if the thermometer is stored or stored outdoors.
- The thermometer should only be used indoors.
- The thermometer is waterproof and can be easily cleaned and disinfected. Clean it before and after each measurement is taken with a cloth dipped in a mild, alcohol-based solution. You can disinfect the thermometer with a suitable solution. The thermometer must never be boiled.
- Keep the thermometer in its case when not in use.
- Remove the batteries if the thermometer is not to be re-used for 3 months.
- LANAFORM cannot be held responsible for any failure of the thermometer if the care and maintenance instructions are not followed.

**FEATURES**

- Rectal, oral and under-arm measurements
- Fast, 10-second measurement process
- Flexible tip
- Audible alarm
- A short beep is heard once the thermometer is started.
- One long beep is heard 10 times when the temperature reading is below 37.8°C (98.6°F).
- 3 short beeps are heard 10 times when the temperature reading is above or equal to 37.8°C (98.6°F).
- Switches off automatically: the thermometer switches off automatically after 10 seconds of inactivity.

**INFORMATION ABOUT THIS THERMOMETER**

The value of the temperature being measured varies depending on the location on the body where it is taken. For people in good health, it may vary between 0.2°C and 1°C (0.4°F and 1.8°F) on various parts of the body.

**Normal temperature ranges between different thermometers**

Temperature values	Thermometer used
Forehead temperature 35.8°C to 37.8°C (96.4°F to 98.6°F)	Forehead thermometer
Ear temperature 35.8°C to 38°C (96.4°F to 100.4°F)	Ear thermometer
Mouth temperature 35.5°C to 37.5°C (95.9°F to 98.6°F)	Ordinary thermometer
Rectal temperature 36.6°C to 38°C (96.8°F to 100.4°F)	Ordinary thermometer

**Normal temperature depending on age**

Age	Temp. °C
0-2 years	36.4 to 38
3 to 10 years	36.1 to 37.8
11-65 years	35.9 to 37.6
>65 years	35.8 to 37.5

**LANAFORM's advice**

Never compare temperatures taken with different thermometers.

Tell your doctor which thermometer you have used to take your body temperature and where on the body. Please also bear this in mind when self diagnosing.

The thermometer may overheat if it is held too long in the hand before taking your temperature. This may skew the result of the operation.

**DESCRIPTION OF THE DEVICE THERMOMETER (FIGURE 1)**

- Display screen
- STOP/START button
- Battery compartment cover
- Probe

**Display (FIGURE 2)**

- Memory button function
- Temperature reading
- Temperature unit
- Low battery indicator

**Installing the batteries (FIGURE 3)**

- Remove the battery compartment cover.
- Insert 1 LR41 battery (supplied) into the compartment, making sure it is facing the right way. The positive (+) and negative (-) polarity symbols are displayed inside the compartment cover.
- Replace the cover.

**Replacing used batteries**

If the low-battery symbol lights up when the thermometer is switched on, replace the battery with a new one.

**SETTING Selection of the unit**

If the thermometer is switched off, switch it on by pressing it continuously for 3 seconds to select a unit of measurement. After 3 seconds, the screen will display either 'C', or 'F'. Start the operation if the unit displayed on the screen does not meet your expectations.

**INSTRUCTIONS FOR USE Measurement process Measuring body temperature**

Remove the protective film covering the battery compartment before using your thermometer for the first time.

- To start the thermometer, press the STOP/START button briefly. A short audible signal indicates that the device is switched on and all the existing icons will appear on the screen for a second.
- To 'C' (°C) and the value of the last measurement taken (memory function) appear one after the other and briefly on the screen.
- The 'C' (°C) measuring symbol then flashes and the thermometer is therefore ready for use.
- The 'C' (°C) measuring symbol then flashes and the thermometer is therefore ready for use. The measuring process ends when the temperature stabilises. An audible signal (10 long beeps) is then heard, the 'C' (°C) constantly flashes when the value of the temperature taken is permanently displayed.

**Storing readings**

The thermometer automatically stores the last reading and it is briefly displayed on the screen before every temperature is taken.

**DECLARATION OF CONFORMITY Declaration of compliance CE 0197**

**Medical device Category IIa.**

We, the company LANAFORM SA, declare that the medical device Digital Thermometer DT-100 LA090113 complies with European Directive 93/42/EEC concerning medical devices and its most recent amendment 2007/47/EC.

**CONTACT INFORMATION The device is manufactured by:**  
Hangzhou Universal / Electronic Co., LTD / 38 Yangjiatang, Sandun / 310030 / Westlake District / Hangzhou China

**The device is distributed by:**  
LANAFORM SA / Rue de la Légende 55 / 4141 Belgium

**Authorised European representative:**  
Name: Prolinx GmbH / Add: Brehmstr. 56, 40239 Dusseldorf, Germany / Tel: 0049-2131 40519680

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**ADVICE REGARDING THE DISPOSAL OF WASTE**

All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.

When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law. Prior to disposal, please remove the battery and dispose of it at a collection point so that it can be recycled.

**Under no circumstances should spent batteries be disposed of with household waste!**

**LIMITED WARRANTY**

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

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Hangzhou Universal / Electronic Co., LTD / 38 Yangjatang, Sandun / 310030 Westlake District / Hangzhou China
Ime: Prolinx GmbH / Add: Brehmstr. 56, 40239 Dusseldorf, Germany / Tel.: 0049-2131 40519680